

GENESIS™

BLOOD COLLECTION MONITOR

MODEL# CM350
CAT# 425-CM350



Keep a close watch on your donation – automatically!


The GENESIS™ CM350 Blood Collection Monitor is a compact, digital whole blood collection scale, ideal for therapeutic phlebotomy and allogeneic blood collection.

Features:

- Electronically monitors and measures blood flow and volume
- Automatically clamps tubing when target volume is achieved
- Rechargeable battery and AC power adaptor
- Automatic calibration
- Audible and visual alarms
- Pole stand sold separately

GENESIS™ BLOOD COLLECTION MONITOR

Designed to support safe collections while maintaining precise control from start to finish.

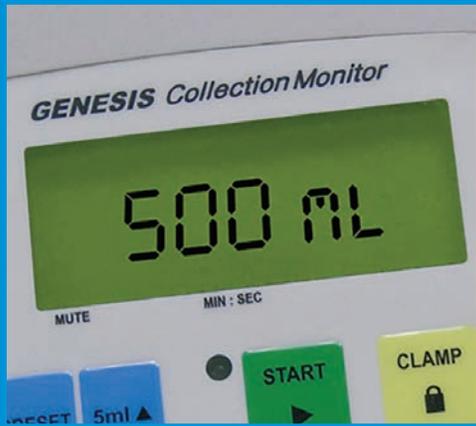

LCD Display provides critical collection information.



Automatically clamps tubing at end of collection.

Ordering Information	Blood Collection Monitor
Model#	CM350
CAT#	425-CM350

Specifications	
Dimensions	6.4"L x 2.4"W x 4.8"H (162L x 60W x 120Hmm)
Weight Display	0 -1199gm Range
Target Volume Range	100 - 600mL Range
Elapsed Time Display	0:00 -19:59 (min:sec)
Power Requirements	Rechargeable NiMH, 7.2V/1.8ah AC Current
Warranty	1 Year

Accuracy	
Weight Tolerance	0 to +3 Grams
Weight Accuracy	+/- 1 Gram
Time Display Accuracy	1 Second

Weight	
Collection monitor	1.50 lb
Collection monitor with pole mount and AC adaptor	2.25 lb

Included Components:

- Aluminum Carrying Case
- Pole Mount
- Rear Hook
- 500g Calibration Weight
- AC Adapter and Power Cord

Optional Accessories:

- Stand (Pole and Pedestal)
- Dry Collection Sets Available

